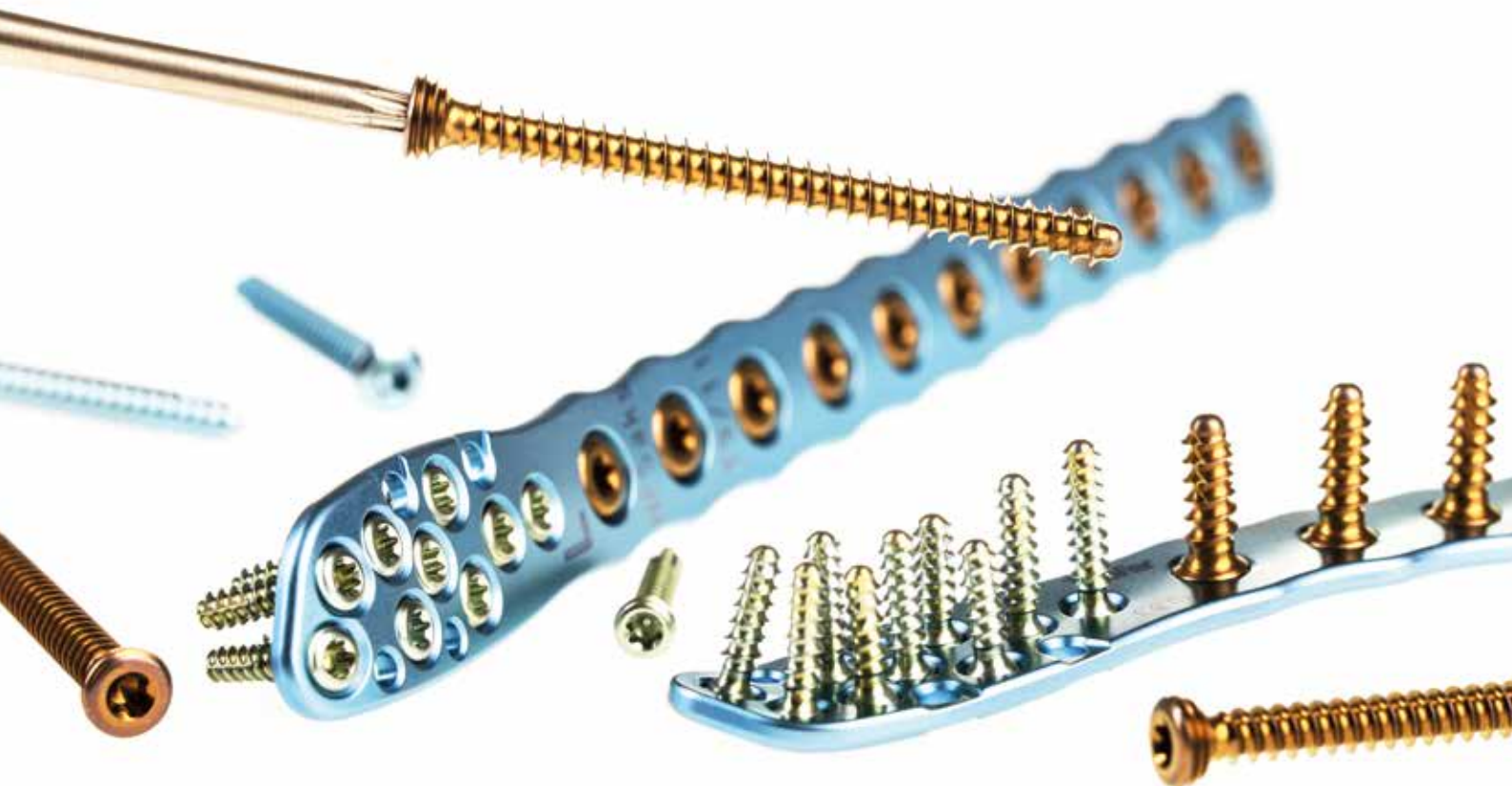
Operationsanleitung  
*Surgical Technique*

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**Warning**

*This description is not sufficient for immediate use of implants and instruments.  
An instruction into handling these instruments by a surgeon experienced with them is strongly recommended.*

|   |              |
|---|--------------|
| <b>Vorwort</b><br><i>Preface</i>  | <b>FLL-1</b> |
| <b>Einleitung</b><br><i>Introduction</i>  | <b>FLL-2</b> |
| <b>Implantatspezifische Informationen</b><br><i>Implant Specific Information</i>                                | <b>FLL-3</b> |
| Indikationen<br><i>Indications</i>  | FLL-3.1-4    |
| Patientenlagerung und Zugänge<br><i>Position of Patient and Approaches</i>                                      | FLL-3.1-4    |
| Plattenspezifische Details<br><i>Plate Specific Details</i>   | FLL-3.1-4    |
| Kontraindikationen<br><i>Contra Indications</i>   | FLL-4.1      |
| Empfehlungen<br><i>Recommendations</i>  | FLL-4.2      |
| <b>Operationstechnik</b><br><i>Surgical Technique</i>   | <b>FLL-5</b> |
| Einbringen des Implantats<br><i>Insertion of Implant</i>  | FLL-5        |
| Temporäre Fixierung<br><i>Temporary Fixation</i>  | FLL-5        |
| Ausrichtung des Implantats<br><i>Orientation of Implant</i>   | FLL-5        |
| Schraubenplatzierung<br><i>Screw Placement</i>  | FLL-5        |
| Entfernung der temporären Fixierung oder der Instrumente<br><i>Removal of Temporary Fixation or Instruments</i> | FLL-5        |
| Wundverschluss<br><i>Wound Closure</i>  | FLL-5        |
| Postoperative Behandlung<br><i>Post-operative Treatment</i>   | FLL-5        |
| Materialentfernung<br><i>Material Removal</i>   | FLL-5        |
| <b>Nachwort</b><br><i>Epilogue</i>  | <b>FLL-6</b> |

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*This document provides information about the handling of HOFER implants and instruments.*

*This operation manual shall be considered as an addition and under no circumstances as a substitute to existing literature about surgical methods within orthopaedics and traumatology.*

*The content shall be regarded as a recommendation for a standardized procedure of how to apply the products without addressing the issues of any further necessary tasks, additional operative actions and possible extensions of the surgical technique.*

*The actual selection of the most suitable implant and its implantation method has to happen exclusively by the surgeon based on his education and the individual diagnostic findings.*

*All illustrations printed here have a purely symbolic character to support the description of the surgical technique and can vary.*

*These operation instructions don't contain any details on the use of the instruments. Corresponding documents are available in the form of*

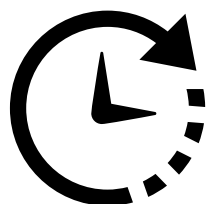
- *Instruction manual for instruments: intra and post-operative handling*
- *Instruction manual for implants (enclosed to each implant)*

*Please note that it is the surgeon's function to identify and characterize the respective injury and its subsequent treatment.*

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*For a safe handling and for the various surgical techniques Hofer-medical is pleased to offer detailed training.*

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Tel.: +43 (0)3382 53388

Mail: [office@hofer-medical.com](mailto:office@hofer-medical.com)

**Please contact our 24/7 service hotline:**

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The **INTEOS**<sup>®</sup> system of the company Hofer GmbH & Co KG (HOFER) is a multi directional and angle stable implant system for small fragments based on the internal fixateur principle.

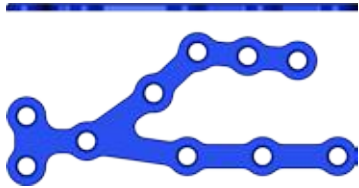
The **INTEOS**<sup>®</sup> system is designed to treat diaphyseal, metaphyseal as well as intra-articular (epiphyseal) fractures of small bones.

All Hofer products result out of a joint development of experienced clinicians and our engineers. This successful cooperation results in providing products to meet the anatomical and functional requirements of the respective sites due to their pre-contoured and low-profile design and to provide an almost unrestricted operative treatment ranging from simple to comminuted fractures.

**General advantages of the system:**

- *Minimization of irritations of ligaments and soft tissue due to the flat plate and screw design*
- *Screw heads countersunk to a large extent*
- *High angular variability of up to 70° with full angular stability*
- *High flexibility of treatment due to sieve design of respective plates*
- *Angular stability in the shaft area: high security in case of bone instability*
- *System added value: compatibility of screws and instruments with many other **INTEOS**<sup>®</sup> plates*
- *3-point support (transformation of bending stress into tensile and compressive stress) is possible*

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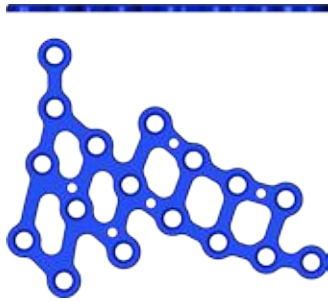
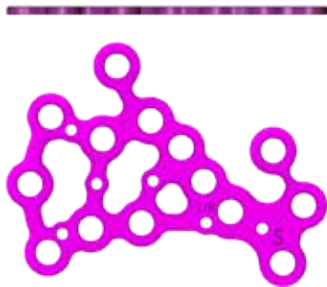
Position of Patient and Approaches

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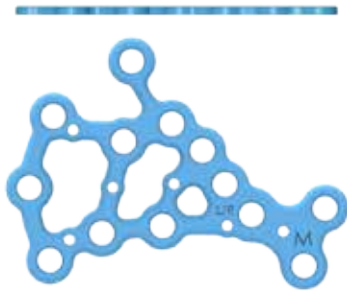
##### Position of Patient and Approaches

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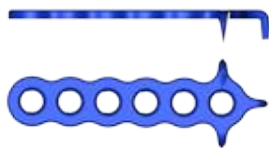
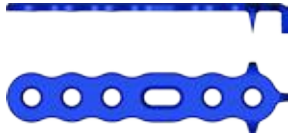
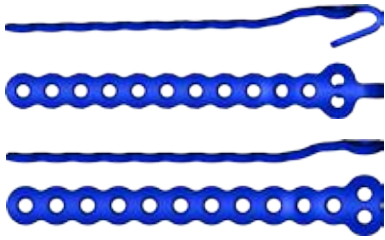
Compatible to INTEOS® 3.0 bone screws  
Recommended screw types:

Bohrergrößen (Standard) für  
Kortikalisschrauben Ø 3,0 2,5 mm  
Drill bit sizes (standard) for  
Cortical screws Ø 3.0 2.5 mm

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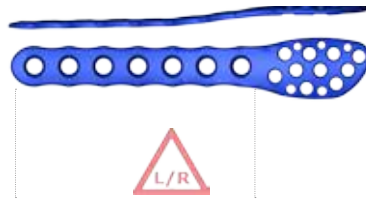
Proximal area  
Compatible to INTEOS® 3.5 bone screws  
Distal area  
Compatible to INTEOS® 2.5 & 3.0 bone screws  
Recommended screw types:

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Drill bit sizes (standard) for  
Locking screws Ø 2.5 2.0 mm  
Cortical screws Ø 3.0 2.5 mm  
Cortical screws Ø 3.5 2.7 mm

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#### **4 - Using the INTEOS® System**

##### **IMPORTANT:**

**Preparatory measures for the use of HOFER implants require a preparation as thorough as possible of the operation field. Nearby nerve fibers and blood vessels require a special caution.**

**Primary an adequate reduction of the anatomical structure must have been carried out before HOFER implants are used.**

##### **4.1 - Insertion of Implant**

After reducing the fractured zone the plate can be positioned. For this step no special insertion device is required. If necessary, a re-shaping of the plate can be performed.

For more details on bending plates please refer to:

- Instruction manual for handling the bending pliers

##### **4.2 - Temporary Fixation**

To temporarily fixate the plate, depending on the plate type, Kirschner wire holes or a slotted hole are available. For the former case the Kirschner wires can additionally be used in a "joy-stick" like fashion, if a further reduction of the respective fragments is required. For the latter case it is recommended to use a not angle stable screw. For information on inserting a screw please refer to the paragraph screw placement (4.4).

##### **4.3 - Orientation of Implant**

The orientation of the implant is carried out according to the present anatomy.

In the case of using the slotted hole: While the screw is not completely tightened, the plate can still be moved to obtain the final position for fixation.

##### **4.4 - Screw Placement**

For placing angle stable or not angle stable screws the technique is the very same. The number of screws, their insertion site and direction has to happen based on the current situation.

For information on drilling pilot holes for the screws please refer to:

- Instruction manual for handling the double drill guide and the drill guide

For measuring the pilot hole depth and for determining the required screw length a depth gauge is to be used.

For more details on hole depth determination please refer to:

- Instruction manual for handling the depth gauge

For more details on inserting screws please refer to:

- Instruction manual for handling the screwdriver

A bicortical screw placement is recommended.

##### **4.5 - Removal of Temporary Fixation or Instruments**

Did the temporary fixation happen using the slotted hole the already placed screw can be used to further fixate the plate to the bone. Therefore, the screw has to be tightened.

##### **4.6 - Wound Closure**

##### **4.7 - Post-operative Treatment**

An intra-operatively applied bandage can support the post-operative analgesia. The starting point for the post-operative functional treatment has to depend on the fracture type and the intra-operatively achieved stability.

##### **4.8 - Material Removal**

In case of secure bony consolidation the surgeon in charge has to decide on the material removal from case to case. For the removal special attention has to be paid to a thorough exposure of the screw heads!

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© 2012 Hofer GmbH & Co KG. All rights reserved.

*Operation Manuals, Handbooks, Handouts and Software are proprietary. Copies, duplicates, translations or conversion of any form in a whole or only in parts are not allowed without a prior written approval by HOFER.*

*Any further rights for the software are specified in the provided licence regulation.*

*The HOFER logo is a registered trademark of the Hofer GmbH & Co KG in Austria and other countries. The rights at other brand and product names mentioned in this document lie with its owners and are recognized hereby. The mentioning of products which aren't from HOFER serves exclusively for information purposes and doesn't represent advertisement of any form. HOFER doesn't take on any liability with regard to the choice, performance or usability of the products.*

*The methods specified in this document, modes of operation and the like were carefully checked. If faults or suggestions for improvement should stand out, however, despite this examination, then we ask you to inform us about these.*

*The operation method represented here introduces merely a possible method for the indication to be treated. It remains to the decision of the respective surgeon to hold on to the described surgical technique or to vary it according to the requirements. Direct and indirect damages as well as consequential damages are therefore excluded in every case.*

*Please do not hesitate to contact us for further information.*

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Ihr HOFER-medical solutions Team

*For us, a partnership is the result of a long and happy relationship in all areas of our work . Reliable bone healing for patients, optimal handling of the systems for surgeons and surgical personnel, as well as the simple preparation of the instruments, constantly strengthen this partnership.*

*For us, "creating" means more than just finding solutions in the form of innovative products.*

*It is the result of high standards, constant development, innovative products and excellent service for patients, surgeons and surgical personnel around the clock.*

*Please do not hesitate to contact us if you have questions about our company, our employees or our production methods.*

*Your HOFER medical solutions Team*