


# LINEAR STAPLER

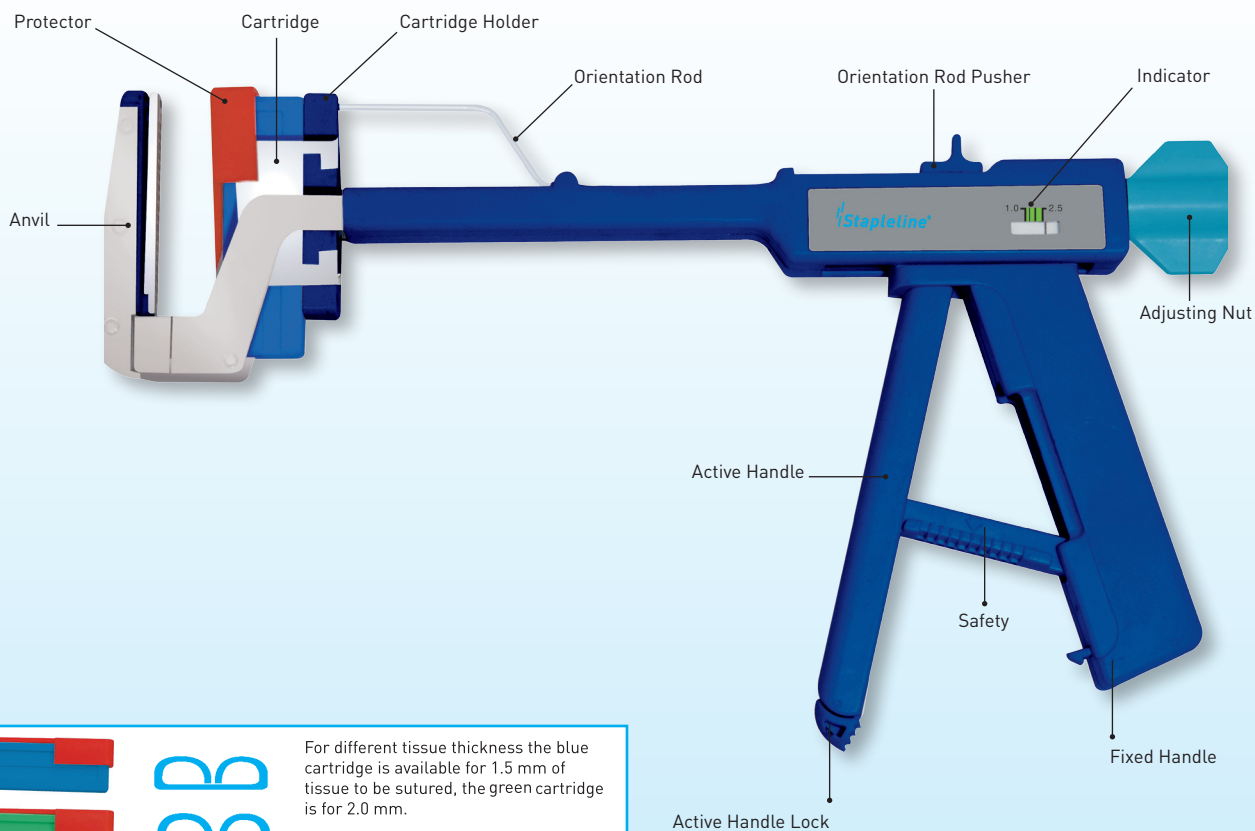

For different tissue thickness the blue cartridge is available for 1.5 mm of tissue to be sutured, the green cartridge is for 2.0 mm.

Safe Suturing by Titanium staples.

The safety can only be released in the right range of tissue thickness.

The adjusting nut enables a comfortable twisting.

Specification of Staples

Wide distance between anvil and cartridge facilitates the insertion of the tissue into the device.

During the surgery, the orientation rod is pushed to ensure the formation of the triggered staples.

Indicator shows clearly the suturing range.

Article	Instrument Order No	Cartridge Order No.	Tissue Thickness*	Leg Height*	Number of Staples	Suture Length*	Cartridge Colour	Sales Unit
FH-Y-30	FH-Y-30-3.8	CLB-30	1.5	3.8	11	32	blue	3 pcs
FH-Y-45	FH-Y-45-3.8	CLB-45	1.5	3.8	15	43	blue	3 pcs
FH-Y-60	FH-Y-60-3.8	CLB-60	1.5	3.8	19	53	blue	3 pcs
FH-Y-90	FH-Y-90-3.8	CLB-90	1.5	3.8	33	92	blue	3 pcs
FH-Y-30	FH-Y-30-4.5	CLG-30	2.0	4.5	11	32	green	3 pcs
FH-Y-45	FH-Y-45-4.5	CLG-45	2.0	4.5	15	43	green	3 pcs
FH-Y-60	FH-Y-60-4.5	CLG-60	2.0	4.5	19	53	green	3 pcs
FH-Y-90	FH-Y-90-4.5	CLG-90	2.0	4.5	33	92	green	3 pcs

\*mm

