

# DISPOSABLE TROCAR



## The trocar is an important instrument in endoscopic surgery.

The Stapleline trocar set consists of 2 sizes of trocar sleeves and obturators, together with adapter sleeves, Veress needle and endobag.

Endobags are available in two sizes: medium and large.

Three sizes of trocars are available: 5 / 10 / 12 mm.

Its properties are as follows:

1. The Stapleline trocar is a sterile single-use instrument.

2. The 10 mm trocar has a spring-mounted protective cover which protects the scalpel.

3. The trocar is marked so that the surgeon can judge the position of perforation protection. During the puncture process, the blade of the scalpel only protrudes 2.5 mm: this prevents any injuries to inner organs.


4. The Stapleline trocar sleeve is gastight and has a CO2 insufflation port.

5. The Stapleline trocar sleeve has grooves which hold the trocar in the tissue, at the same time they stabilize the sleeve which thus has a smaller diameter. In this way, the wound can be kept as small as possible. There are 2 lateral openings to fasten the trocar sleeve during surgery.



Trocar Set 1	Order number
2 x Trocar sleeve 5 mm	C04-02-Z05B
2 x Trocar sleeve 10 mm	
1 x Trocar insert 5 mm	
1 x Trocar insert 10 mm	
1 x Reducer	
1 x Veress needle 118 mm	
1 x Endobag 10 x 13,5 cm	

Trocar Set 2	Order number
2 x Trocar sleeve 5 mm	C04-02-Z06B
1 x Trocar sleeve 10 mm	
1 x Trocar sleeve 12 mm	
1 x Trocar insert 5 mm	
1 x Trocar insert 10 mm	
1 x Trocar insert 12 mm	
1 x Reducer	
1 x Veress needle 118 mm	
1 x Endobag 10 x 13,5 cm	

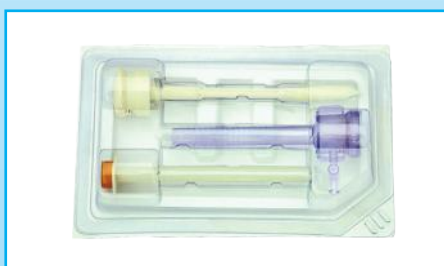

Trocars singly	Order number
10 x 5 mm (w/o.S.)	C04-02-05-10
10 x 5 mm (w.S.)	C04-02-05-10-2
10 x 10 mm(w.S.)	C04-02-10-10
10 x 10 mm (w.S.) long	C04-02-10-10-L
10 x 12 mm (w/o.S.)	C04-02-12-10
10 x 12 mm (w.S.)	C04-02-12-10-1
10 x 12 mm (w.S.) long	C04-02-12-10-L

Separate Items	Order number
10 x Reducer	B04-10-1004
10 x Veress needles 118 mm	B04-10-1005
10 x Veress needles 150 mm	B04-10-1007

(w/o.S.) = without Scalpel / (w.S.) = with Scalpel

